



## New Spectra Iris

Adjustable aperture for all pupil sizes

Keeler



## New Spectra Iris

### The little indirect now has even bigger benefits!






Spectra Iris has been specifically designed as a compact, lightweight and portable indirect with big benefits for professionals on the move. It can be tucked away in its neat carrying case or hung around your neck when not in use.

Plus now, just like the iris in your eye, you can adjust the size of the aperture at the touch of a slider.

### LED long lasting illumination

No more bulbs just long lasting white homogenous light patch that can be used for up to 4 hours continuous on time.




-  Adjustable aperture for all pupil sizes
-  Portable and convenient
-  LED long lasting white illumination
-  Lightweight and comfortable
-  Stylish, comfortable, attractive frame

### **NEW** Adjustable Aperture Slider for all pupil sizes




Now you can adjust the aperture simply by moving the slider.

Slide to the right to increase aperture size and to the left to decrease. The aperture range is from 20mm to 60mm.

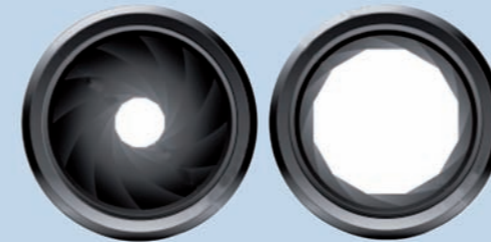

*New Adjustable aperture for all pupil sizes from 20mm to 60mm (not to scale)*

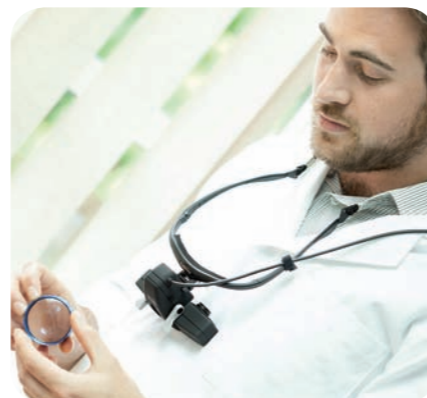

### Keeler Sport Frame

Spectra Iris comes on the new Keeler Sport Frame. This wrap around frame is not only attractive but provides maximum comfort and balance.

### Flip up optics

The whole optical box and light pod can be flipped up allowing direct eye contact when talking to patients or for writing notes.



### Filters and apertures

P D Range: 48 - 76mm



- 1 Aperture slider
- 2 Filter selection control
- 3 Mirror angle control
- 4 Keeler Sport Frame
- 5 Flip up optics

## Variable PD

No restrictions. Just adjust your PD when you need it. It has a wide range between 48 to 76mm.

## Great Lightweight Lithium ion Battery

The small compact Lithium ion battery simply clips on to your belt and can be stored in its charger when not in use.

It provides up to 4 hours continuous use on time.




## Carrying case

A little case for a little indirect.



## Fits neatly and simply over your spectacles

Spectra Iris can fit over the top of your spectacles so no need to remove your correction.

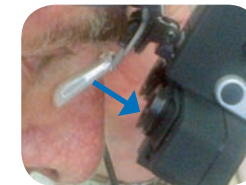

Old hinge



New hinge

## New single movement hinge

The new hinge allows the Spectra Iris to sit closer to your eyes for a better view and a more comfortable fit.

## Part Numbers

Spectra Iris Frame Size	Part Number
Small	1205-P-2012
Regular	1205-P-2010
Large	1205-P-2011

Optional Accessories	Part Number
Spectra Iris Carrying Case	3412-P-7002
Spectra Iris Teaching Mirror	1205-P-7000
Spectra Iris Battery	1919-P-5215
Spectra Iris Lithium Charger	1941-P-5385

Spectra Iris systems are supplied with a lithium battery, charger, indirect and Black Keeler Sport Frame.



### Keeler Ltd.

Clewer Hill Road, Windsor  
Berkshire, SL4 4AA

Tel: +44 1753 857177

Fax: +44 1753 827145

Email: [info@keeler.co.uk](mailto:info@keeler.co.uk)

Website: [www.keeler.co.uk](http://www.keeler.co.uk)

### Keeler Instruments Inc.

456 Parkway, Broomall  
PA 19008 USA

Tel: 1 610 353 4350

Fax: 1 610 353 7814

Email: [keeler@keelerusa.com](mailto:keeler@keelerusa.com)

Website: [www.keelerusa.com](http://www.keelerusa.com)

# Keeler