

Keeler
3222 Phoenixville Pike, Bldg. 50
Malvern, PA 19355
800.523.5620 | 610.353.4350


Extra Large Slit Lamp Breath Shield

Ref: EP39-80753



With the current global health crisis, there has been a huge demand from healthcare professionals for products to increase their protection from airborne pathogens. As such Keeler is introducing an "Extra Large Breath Shield", that will introduce extra protection between the clinician and patient when using a slit lamp.

The breath shield will fit any Keeler slit lamp, being 35cm x 30cm (13.7" x 11.8") gives great protection while still allowing easy access with a Volk lens. The shield is made of 1.5mm clear plastic - making it light but also sturdy.

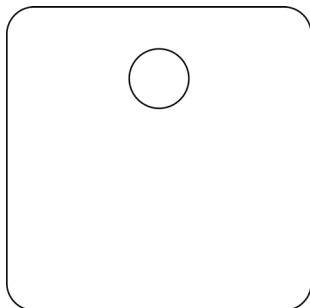


Extra Large Breath Shield

Part Number: EP39-80753

Standard Breath Shield

Part Number: EP39-80242

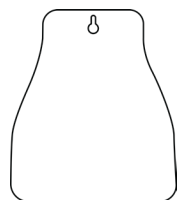

Extra Large Breath Shield Fitting

The breath shield is fitted by removing the binoculars (with the securing knob), and placing the shield over the dovetail joint.

The shield should fit any slit lamp where the outside diameter of the dovetail joint is less than 70mm.

